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


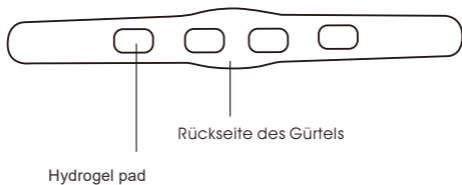
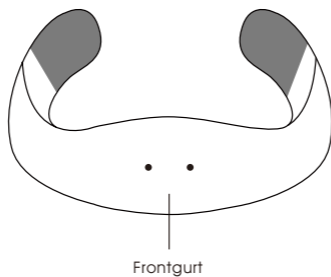
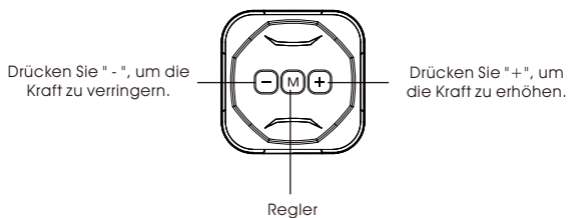
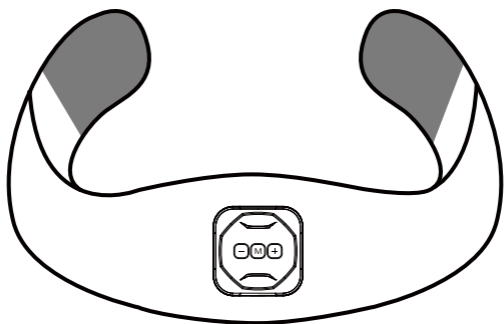
Rechargeable Style Slimming Massage Belt

User Manual

## 1. Product Specification

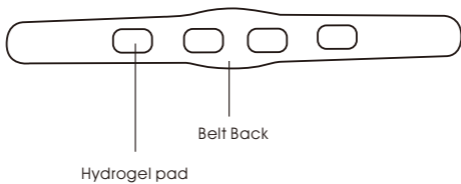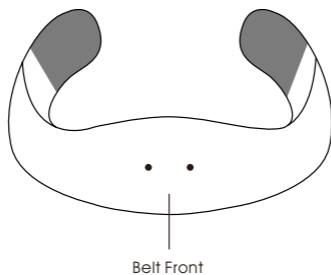
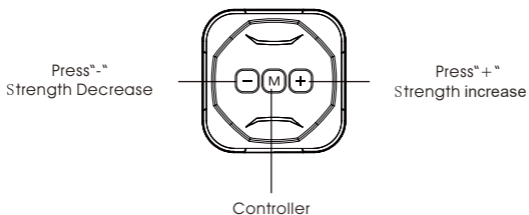
Product Name :	Rechargeable Style Slimming Massage Belt
Model:	FIT57
Working voltage:	DC3.3V
Working Environment temperature:	10-40°C
Charge Voltage:	DC5V
Battery Type:	Polymer lithium battery
Package size:	228 x 210 x 46mm
Executive standard :	GB4706.1-2005 GB4706.10-2008
Accessories:	Controller x 1 ;User manual x 1 ;Gel pad x 1 ; USB cable x 1

## 2.How to Use

1. Pls fix the controller to the belt.

2. Place the belt into the right location of your abdomen

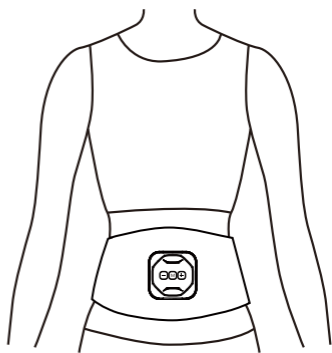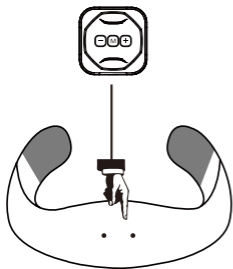
3. Press M Button to turn on the controller, the indicate light will be turn on. It will change different mode when press M button, total have 6 modes.

4. Press + key to start the work, and then press + key to adjust the intensity. Press + one time the intensity will increase one grade. Press - key one time, it will decrease one grade.

5. When you want to turn off, please press - key for several times until the power turn off.

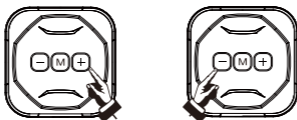Mode 1  
One white light lights up

Mode 2  
Another white light lights up

Mode 3  
Two white lights light up

Mode 4  
One blue light lights up

Mode 5  
Another blue light lights up


Mode 6  
Two blue lights light up


## Operation and maintenance of the Hydrogels Pads

### Main points:


If the adhesion method is not correct the adhesion force with the skin may be decreased and the function cannot be performed normally. Or shorten the shelf life due to excessive consumption.

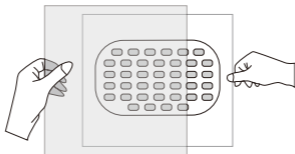1. Take the bottom plate from the bag.
2. Fix the controller to the gel pad.
  - Install the gel pad to the correct location of the controller.
  - Pay attention to the pad for prevent the two sides to glue it together when install, because the gel pad is very soft.

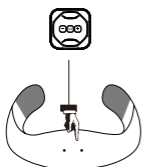
Belt Front

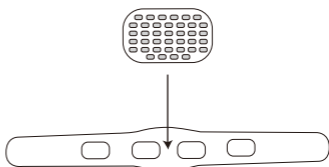
Belt Back

3. Peel off the gel pad from the bottom plate.

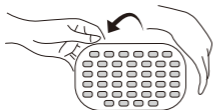4. Put the pad to the use area when you installed well the products and use the pad into the body for exercise.

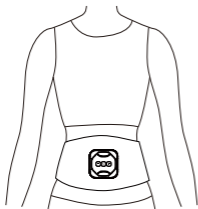Please follow the following steps for maintenance to avoid the adhesion surface of the gel pad difficult to attach due to dirty.

1. Please dip your finger into the water droplets and wipe the adhesive surface of the hydrogel to clear the dirty part.
2. After fully draining the water, place the adhesive surface (adhesive facing up) in a well-ventilated cool place to dry.
3. Dry naturally until there is no moisture on the surface, then paste the hydrogel back to the bottom plate

Important: Don't take long or frequent maintenance, or it may damage the adhesive surface and cause shorten the use lifetime.

## Safety Warning

Not to be used with equipment in the following conditions as it may lead to accidents or malfunctions to the machine, and /or serious bodily damage.

- \* The use of heart regulator and Internal transplant electronic Instruments for medical use.
- \* The use of maintain life equipment such as artificial heart and lung for medical use.
- \*The use of Electrocardiogram and other wearable for medical use.


Not to be used with person in the following conditions as It may lead to accidents or skins and other bodily damage.

- \*Patients who had heart disease or high blood pressure.
- \*Pregnant women and postpartum women.

Please consult your doctor before use If you had following conditions as It may lead to accidents or skins and other bodily damage.

- \*The area of Irritation on the skin.
- \*Patients with blood circulation disorders.
- \*Patients who taking the medicine.
- \*Patients with malignant tumor.
- \*Patients who had surgery recently.
- \*Intuitive disorder caused by high terminal circulatory disorder, such as diabetes
- \*Patients with fever 38 degrees above, patients need take good test.
- \*Acute diseases such as spinal fracture, sprain, muscle strain, etc.
- \*Patients with abnormal osteoporosis.
- \*Patients with abnormal blood pressure.
- \*Patients with abnormal body.
- \*Patients being treated at a medical facility.

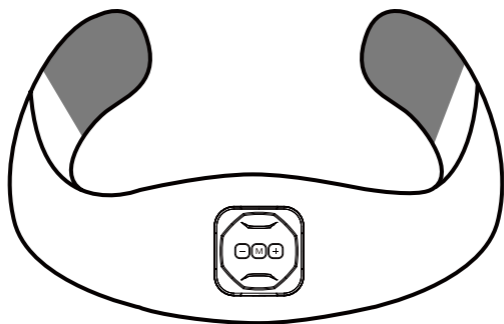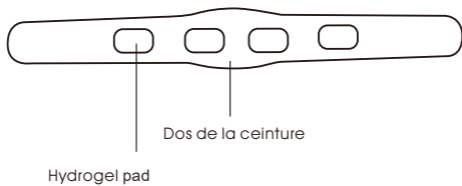
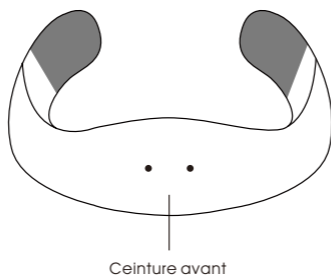
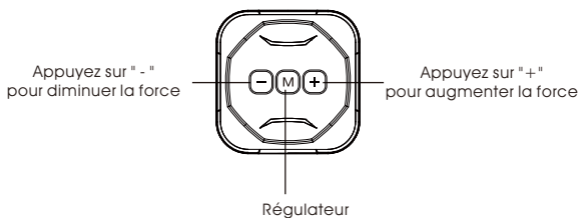
Equipment and remotes should not place or near bathroom or close to the fire as It may cause electric shock or malfunction to the machine.

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